edel
white

SONIC GENERATION®

WINNER

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Instruction manual

Manuale d'uso

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HERZLICHEN GLÜCKWUNSCH

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For safety reasons

Warning:

- If you have had oral or gum surgery lately, please consult your dentist. Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues.
- Perform and follow through breathing exercises responsibly; it's under your control. Always do breathing exercises in a sitting position, for your own protection. Stop breathing exercises as soon as you feel pain, discomfort or light-headedness. Never force yourself to the point of discomfort. If you take medication or are under medical treatment, always consult your physician before practicing any breathing exercises.
- Do not take the device apart, (Exception: Disposal).

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe manner, and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the appliance.


Caution:

- The supply cord cannot be replaced. Dispose of the device if the cord is damaged.
- Keep charger away from water or other liquids! Do not leave it near a bathtub or a basin.
- Keep charger away from hot surfaces and direct sunlight if connected to the socket.
- The Sonic Generation® WINNER toothbrush complies with safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device.

Please read this instruction manual carefully before using the device and save it for future reference.

THE TECHNOLOGY OF SONIC GENERATION® WINNER

Sonic Generation® WINNER is a rechargeable hydro-active toothbrush with sonic-technology for a 100% all- round tooth cleaning. Surges of the sonic technology create pulsating flows of fluids, which easily penetrate even the narrowest dental spaces. It eliminates debris and discolorations all around every single tooth for naturally white and healthy teeth.

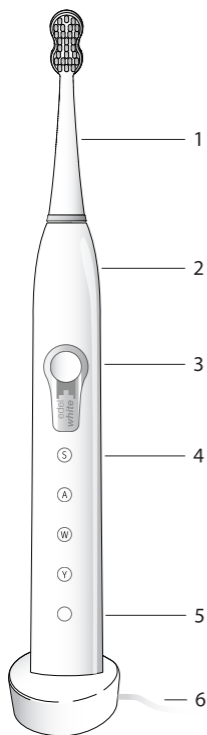
- 1.) Dual Clean brush head with travel cap
- 2.) Handle with brushing timer and quadrant timer
- 3.) On/off button and Mode Selection
Power-on: Press button briefly.
Mode Selection with memory to change brushing mode: Press the button again within 3 seconds after turning on the device. Repeatedly press the button to select the desired function.
Power-off: Press the button at least 3 seconds after turning on the device. Alternatively, press the button repeatedly until the device turns off (after clean mode YOGA).
- Clean mode SENSITIVE (S):**
Upper LED-light glows. Gentle cleaning for sensitive teeth and delicate gums.
- Clean mode ACTIVE DEEP (A):**
Upper two LED-lights glow. Thorough plaque removal for healthy teeth and strong gums.
- Clean mode WHITENING (W):**
Three LED-lights glow. Natural tooth whitening by removing tough stains.
- Clean mode YOGA (Y):**
Four LED-lights glow. Intensive cleaning action and „breathing trainer 4-4-8“ in one.
- TRAVEL LOCK:**
Press for 4 sec.: locks and unlocks the button, and prevents accidental operation during travel.
- 4.) LED display: Indicates the selected clean mode.
Pulsating LED-lights: Indicates the state of charge.
- 5.) LED display of travel lock
- 6.) Induction charger for rechargeable Li-Ion battery, micro-USB cable (1 m) and power adaptor with automatic voltage adjustment 100-240 V plus tray with 3 satellites for brush heads and 3 mounts for interdental brushes

OPERATION

Align the brush head so the bristles face the front of the handle. Firmly press the brush head down on the metal shaft without twisting until it clicks into place.

Use the USB cable and connect the charger with the power adapter. Plug the adapter into a socket and place the hand piece onto the charger. The LED display will light up for about 10 sec. Afterwards all LED-lights will switch off. The Sonic Generation® WINNER toothbrush is charging.

Pull the brush head off the handle without force for replacement. Don't twist the brush head when removing it from the handle as it may damage the motor. Don't immerse the toothbrush handle or the charger under water and do not leave the handle on wet surfaces. Dry the toothbrush handle off with a soft cloth after brushing.

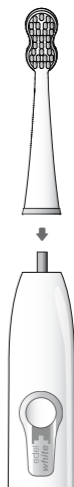CHARGE - LEVEL INDICATOR

The charging takes about 12 hours. Charge the Sonic Generation® WINNER toothbrush completely before first use.

The LED display will blink for 10 sec. if the device is removed from or placed on the charger, or switched off.

Four blinking LED-lights indicate a fully charged battery. If one single light is blinking the battery is charged less than 25%.

Don't keep the toothbrush handle on the charger base all the time. Charge the toothbrush only when the charge-level is low, to maximise battery life.

USAGE OF THE SONIC GENERATION® WINNER

Wet the bristles and apply a small amount of toothpaste.

Gently place the brush against the teeth at a slight angle towards the gum line.

The toothbrush remembers the last chosen mode and will start automatically by pressing the on button in this mode (Memory). Repeatedly pushing the button within 3 seconds after turning the device on will change the cleaning mode:

SENSITIVE - Modus: Upper LED-light glows

ACTIVE DEEP - Modus: Upper two LED-lights glow

WHITENING - Modus: All Three LED-lights glow


YOGA - Modus: Four LED-lights glow

Move the toothbrush in an up-and-down manner from tooth to tooth, without exertion; let the brush do the work. Rest 2-3 sec. on one tooth before moving to the next one. Make sure that the longer bristles penetrate the dental spaces between teeth.

For each quadrant, move the brush first along the internal tooth side over to the chew surface, and at last along the external tooth side.

After 30-40 sec., the toothbrush signals to change the quadrant by a short interruption in its vibration. Repeat this cleaning process for each quadrant. When the time has run out the device automatically deactivates. If needed, it is possible to pursue the brushing by simply reactivating the device.

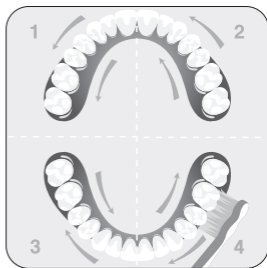Important: do not exert pressure on the toothbrush. Gentle surface contact along the tooth alignments will lead to better cleaning results and maintain the health of your teeth and your gum line.

Your Sonic Generation® WINNER is also safe to use on: braces, implants and dental restorations (fillings, crowns, veneers).

Replace Sonic Generation® WINNER brush heads every 3 months for hygienic reasons. Brush heads wear out sooner when used on braces.

Ask your dentist, dental hygienist or assistant after every appointment if and how you are supposed to use your sonic toothbrush. Dental professionals are always competent for follow-up questions.

CLEAN MODE

Select the different modes by repeatedly pushing the button. The Sonic Generation® WINNER remembers the last chosen mode so that it does not have to be selected the next time you activate the device (**Memory function**).

Clean mode SENSITIVE

By selecting SENSITIVE mode the upper LED-light glows. The sonic waves that move with a speed of up to 42.000 strokes per minute allow the bristle tips to vibrate gently and penetrate even the smallest dental spaces. This extremely smooth and efficient way of cleaning is very gentle on enamel, gums and exposed dental necks. We recommend the SENSITIVE mode especially to people with sensitive teeth, exposed dental necks, periodontitis, implants, bridges, and other dental replacement solutions.

Clean mode ACTIVE DEEP

By selecting ACTIVE DEEP mode the two upper LED-lights glow. The surges of the sonic technology create a pulsating flow of fluid which easily penetrates even the narrowest dental spaces, and neutralizes harmful debris, even in hidden niches. This mode assures ideal cleaning and health of your teeth. For healthy teeth it is recommended as the basic mode.

CLEAN MODE

Clean mode WHITENING

By selecting WHITENING mode three LED-lights glow. Increased brushing movements on the tooth surfaces allow increased cleaning. Teeth are freed of unpleasant discoloration. We recommend the WHITENING mode especially to cases of discoloration as a result of nicotine, wine, tea or coffee, and for occasional extra-thorough teeth cleaning.

Clean mode YOGA

By selecting YOGA mode all four LED-lights glow. SONIC GENERATION® WINNER supports you twice a day with breathing exercises while you clean your teeth. A good breathing rhythm has direct impact on our well-being and our whole bodies. For practice, SONIC GENERATION® WINNER suggests timing your breaths „4-4-8“ (breathe in for 4 seconds - hold your breath for 4 seconds - breathe out for 8 seconds). Repeat as you move the brush around your mouth. One quadrant of your teeth can be brushed and cleaned with each 4-4-8 cycle. Stop brushing while the device pauses for 4 seconds as you hold your breath. The cycle is repeated 12 times on each surface that's cleaned before the device turns off. Please note the warnings in this manual.

Quadrant timer

The quadrant timer is a 30 second interval timer that briefly interrupts the brushing to signal the recommended change to the next quadrant. It reminds you to brush the four-quadrants of your mouth equally and thoroughly, for a better overall-cleaning.

Brushing timer

Dentists recommend a minimum of 2 minutes brushing time. After the recommended brushing time, depending on when the selected mode; has been reached, the Sonic Generation® WINNER deactivates.

TRAVEL LOCK

The Sonic Generation® WINNER- travel lock; locks and unlocks the device in order to prevent accidental operation during travel.

Locking the device: In order to lock the device, push the button and hold it for about 4 seconds, until the LED display of the travel lock (5.) blinks. The device is now safe from accidental activating.

Unlock: To unlock the device push the button for about 4 sec. The LED-lights show by blinking that the device is unlocked. An alternative is to put the toothbrush back on the charger base. This unlocks it automatically.

QUESTIONS AND ANSWERS

I have sensitive teeth and gums. What is the best way to brush with my Sonic Generation® WINNER?

If you have sensitive teeth and gums, plan for a 1-2 week adjustment period to your Sonic Generation® WINNER. During this time, we recommend using the SENSITIVE mode.

How can I best remove discolorations?

We recommend the WHITENING mode in order to rid teeth of discolorations.

What do I have to consider when I brush my teeth with Sonic Generation® WINNER if I am wearing braces or if I have dental restorations (bridges, implants, crowns or veneers)?

In this case use the ACTIVE DEEP mode to remove plaque and bacteria. In case of braces, place the brush gently between your brackets and gums and move the brush in a light circular motion along your teeth to keep both brackets and teeth free from debris.

I have recently had oral surgery.
Can I use Sonic Generation® WINNER?

Always discuss your oral hygiene requirements with your oral surgeon or dentist and follow the given instructions after an operation. As soon as your teeth and gums heal to permit brushing, you can proceed using Sonic Generation® WINNER.

You will get further information at:
www.edel-white.com/sonic

Sonic Generation® WINNER brush heads

edel+white® SOFT DUAL CLEAN

The feather-fine Konex® bristles enter even the smallest dental spaces. Use for sensitive teeth and gum problems.

edel+white® WHITENING DUAL CLEAN

The patented polishing toothbrush with embedded rubber balls erases discolourations and gives them a brilliant white appearance without scouring.

edel+white® SINGLE CLEAN

The feather-fine fronds of this innovative sulcus brush reach further and clean all around brackets, implants and closely spaced teeth.

Cleaning the appliance

Rinse the brush head and bristles after each use with running water. Remove the brush head from the handle and rinse the connection at least once a week. Always disconnect the charger from the socket, wipe it with a cloth and only plug it into the socket if fully dried. If you are not going to use the appliance for an extended period of time, remove the main plug from the wall socket. Then clean the appliance and store it in a cool and dry place.

Warranty

edel+white® guarantees its products for two years after the date of purchase. This warranty is the same in every country in which this device is sold by edel+white® or another authorized dealer.

This warranty does not cover:

Defects due to inappropriate use, regular wearing off or handling, especially of the exchangeable brush, but also flaws that insignificantly affect value or functionality of the device. The warranty runs out if being repaired by unauthorized services or if opening the body housing.

In order to claim warranty, please send the entire device with pro-

of purchase to the official edel+white® partner in your country. If you do not know your official edel+white® partner, please visit our homepage www.edel-white.com/helpline or contact us directly in Switzerland: +41 (0)44 520 0417 (The usual costs of a phone call into a foreign country will be charged).

TECHNICAL SPECIFICATIONS

Mains voltage: 100-240 V AC , 50 Hz, 2.5 W, Class II, IPX7

Battery: rechargeable Li-ion-battery

Operating time: 2 weeks at 2x2 minutes per day

Initial charging time: 12 hours

This device complies with the guidelines of the EU (EMC/RoHS) and low voltage. The power adaptor is GS certified.

Disposal

The battery must be removed from the appliance before it is scrapped. The appliance must be disconnected from the power supply when removing the battery. Dispose of the discharged batteries correctly. Special containers are provided for disposing of rechargeable batteries at shops selling them or public collection points. When disposing the main unit, affix scotch tape or label to the electrodes of the battery. Do not throw the rechargeable battery in fire.

Correct disposal of this product:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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edel+white® and Sonic Generation® are internationally registered trademarks of Scanderra GmbH, Zurich, Switzerland. All rights and modifications reserved.