

## Air Scaler

# TA-200

4HA / 2HA / KA


A single multi-purpose handpiece for scaling and endodontic procedures.

### More efficient and rapid



A versatile handpiece plays double role as air scaler and endodontic instrument by changing attachments. With the vibration by compressed air, TA-200 enables much more efficient and rapid root canal cleaning and enlargement as well as tartar removal than manual treatment.

### Controllable power output

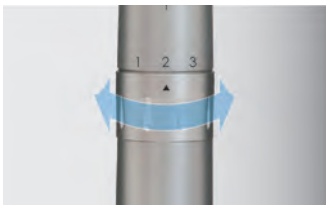

TA-200 is equipped with a power controller. The best power output can be obtained according to the purpose of use, and to the attached instrument.

Scaling

Endodontic



## Specifications

Model	TA-200-4HA	TA-200-2HA	TA-200-KA
Connection	ISO standard 4-Hole	ISO standard 2-Hole	MULTIflex® coupling
Frequency	5,000 to 9,000 Hz		
Air pressure	0.3 MPa ( 44 psi, 3.0 bar, 3.1 kgf / cm <sup>2</sup> )		
Length (without scaling tip)	110 mm (4.3 inch)	104 mm (4.1 inch)	109 mm (4.3 inch)
Weight (without scaling tip)	74 g (2.6 oz.)	70 g (2.5 oz.)	73 g (2.6 oz.)

MULTIflex® is a registered trademark of KaVo Dental GmbH.

Specifications are subject to change without prior notice for improvements.

### Manufacturer : MICRON CORPORATION

2-17-7, Ikegami, Ota-ku, Tokyo, 146-0082, Japan  
 Tel : +81 3 3755 0396, Fax : +81 3 5747 5396  
<http://www.micdent.com/> E-mail : [contact@micdent.com](mailto:contact@micdent.com)

### European Authorized Representative : Advena Ltd.

Tower Business Centre, 2<sup>nd</sup> Flr., Tower Street, Swatar, BKR 4013 Malta