



SMART-CRUISE SMART CHOICE



Zhuhai Siger Medical Equipment Co., Ltd

Add: Building 2, No.1 Chuangxin Yi Road,
Tangjiawan Town, Zhuhai City,
Guangdong Province, P. R. China.

Tel: +86-756-3881630 Fax: +86-756-3881028

Zip Code: 519080 Website: www.siger.cn E-mail: scaler@siger.cn




PROFESSIONAL PERIODONTAL

TREATMENT DEVICES

NP I • NP PRO • NP COMPACT

Perfect Combination of Efficiency and Comfort

- Advanced smart cruise system
- Automatic frequency tracing
- Real-time feedback of load
- On time power compensation

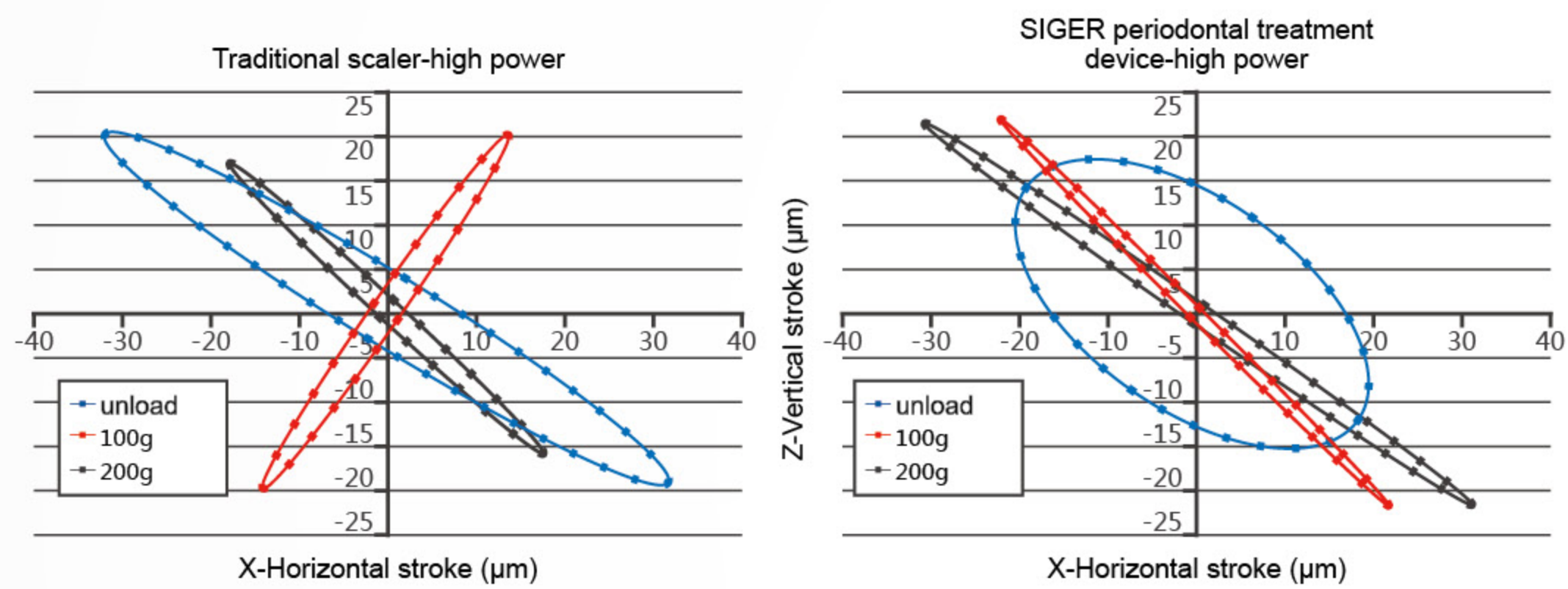

- Elliptical stroke pattern.
- 30 HRC medical titanium tip, low hardness and less injure.
- Frequency: 32KHz
- Vibration Amplitude: 20-90µm
- Treatment depth reaches 14mm.



Before and after treatment

- Less pain treated by Siger periodontal treatment device.



Treated by traditional ultrasonic scaler



Treated by Siger periodontal treatment device

- 600ml Water tank(NP Compact/NP Pro).
- Measurement line on the water tank makes it easy for medical solution matching.
- Hydrogen peroxide, chlorhexidine etc are applicable.



NP COMPACT

Multifunctional Painless Periodontal Treatment Device

SMART-CRUISE, SMART CHOICE

- Advanced smart cruise system.
- Perio and Endo function.
- Two kind of water supply system.
- Tempered glass touch panel.
- Fault indication.
- LED lightening mode optional.
- Autoclavable detachable handpiece.
- Medical titanium tip.

Multifunction touch panel

- Sensitive control, easy cleaning.
- Water proof glass material.
- Perio and Endo function.
- 10 grade power adjustable, high efficiency combined with comfort.
- Cavitation effect and acoustic microstreaming effect.
- Fault indication.
- One-button cleaning.



- LED handpiece, better view during treatment.
- Lightening mode optional.


- Two types of water supply system, medical solution applicable under auto water supply mode.
- Leak-proof water tank design, easily operation.



- Handpiece holder with lock button, anti-maloperation.



Professional perio and endo treatment kit

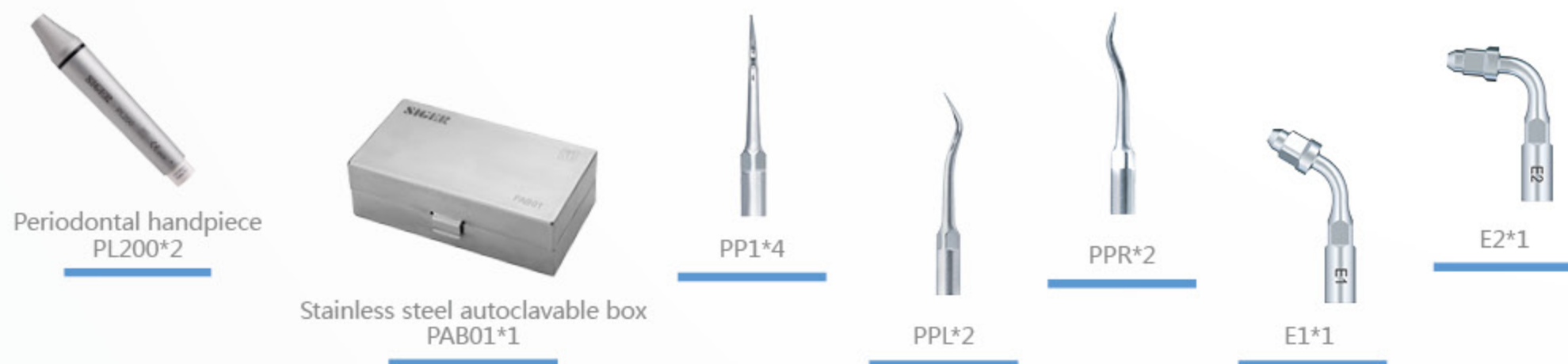

NP PRO

Painless Periodontal Treatment Device

- Advanced smart cruise system.
- Tempered glass touch panel.
- Fault indication.
- Autoclavable detachable handpiece.
- Medical titanium tip.



Professional periodontal treatment kit



Periodontal handpiece
PL200*2



Stainless steel autoclavable box
PAB01*1



PP1*4



PPL*2



PPR*2

SMART-CRUISE, SMART CHOICE


Multifunction touch panel

- Sensitive control, easy cleaning.
- Water proof glass material.
- 10 grade power adjustable, high efficiency combined with comfort.
- Fault indication.



- LED handpiece, better view during treatment.
- Lightening mode optional.

- Medical solution applicable.
- Leak-proof water tank design, easily operation.

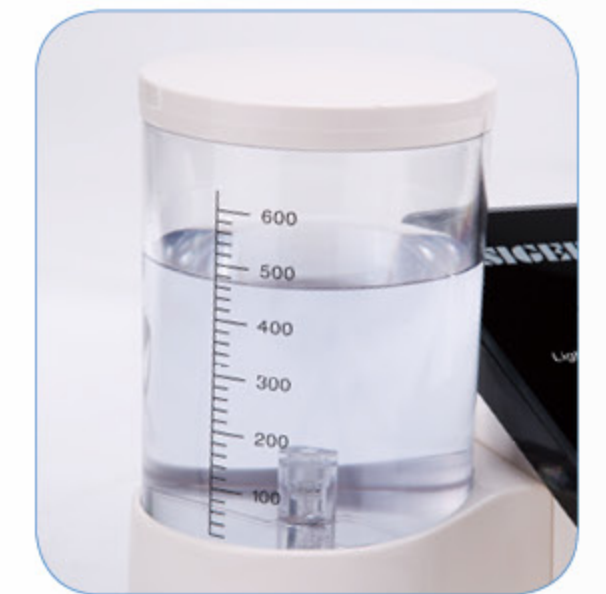

- Handpiece holder with lock button, anti-maloperation.





- Advanced smart cruise system.
- Autoclavable detachable medical handpiece.
- LED handpiece provide better view during operation



- LED handpiece, better view during treatment.
- Scaler tip is with elliptical stroke pattern during operation.
- Medical titanium tip, low hardness and less injure.

Professional periodontal treatment kit



Periodontal handpiece
PL200*1



Stainless steel autoclavable box
PAB01*1



PP1*1



PPL*1



PPR*1

- Tidy and convenient, comfortable and efficient.



Tips application



PP1 (Straight)

- Applicable to the Lingual(lip) side of mandibular(upper) anterior teeth.



PPL (Left curved)

- Applicable to the Lingual(lip) adjacent side of mandibular(upper) anterior teeth.



PPR (Right curved)

- Applicable to the Lingual(lip) adjacent side of mandibular(upper) anterior teeth.



PPB (Ball)

- Applicable to the scaling and planning of root furcation area.



Function description

● Standard / Not available

	NP COMPACT	NP PRO	NPI
Image			
Perio	●	●	●
Endo	●	/	/
Touch panel	●	●	/
LED option	●	●	/
Self-container	●	●	/
Water inlet	●	/	●
Handpiece holder lock button	●	●	/
Fault indication	●	●	/

Spare parts



Icon

